CL-300 Computerized Lensmeter





CL-300 Computerized Lensmeter

With a new UV Measurement Function, the CL-300 provides comprehensive data to both eyeglass wearers and optometrists. Its compact, ergonomic design and intuitive touch screen panel ensures that measurements can be easily and quickly obtained.


Features

- LCD Touch Panel
- Compact Slim Body
- >> UV Transmittance Measurement
- » Green Measurement Light Beam
- Automatic Mono- & Multi-focal Detection

UV Transmittance Measurement Function

The CL-300's new UV Measurement Function provides information on the ultraviolet transmittance in the range of 0% to 100%, providing reliable measurement results to both eyeglass and sunglass wearers.

* UV Transmittance Measurement is performed with the light beam of 365 nm wavelength.

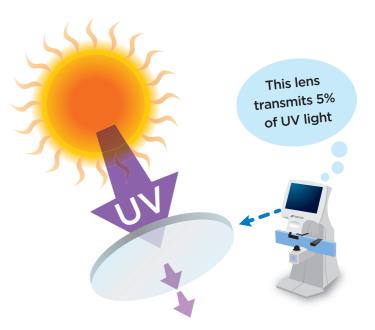

Multi-focal Lens Measurement

Not only single lenses but also multifocal lenses can be easily measured. The CL-300 also features automatic detection for single and progressive lenses.




example: Progressive Lens Measurement

The image shows a near point measurement screen. With simple operation, the ADD value can be obtained by following the guide instructed on this progressive measuring mode.

Green Light Beam

The CL-300 incorporates the green light reading beam, which is also known as e-line. Due to the wavelength of 546 nm, the e-line can be well-sensed by the human's eye and it is determined as a standard by ISO 7944. Using green light beam in our instrument enhances lens measurement precision and reduces the chance of error in the measurements.

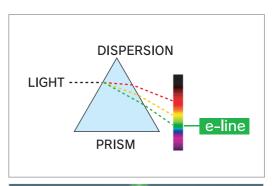



Improved Lens Support

Due to its new design, the lens support of the CL-300 is 6 mm more in height compared to the conventional model. This change enables smooth and trouble-free measurement of high-curve lenses, preventing from troublesome contact of the lens edge with external parts of the lens support.

Large 5.7 inch LCD Monitor with Wider Directional Angle


The wide LCD monitor covers a 160 degree angle of observation allowing clear, compatible readings from virtually any position regardless of the operator's height.

White Ink Cartridge Standard

The newly featured white marking ink provides clear marking even on lenses with anti-staining and anti-fog coating. Due to these favorable qualities, TOPCON made the white marking ink standard in the new CL-300.





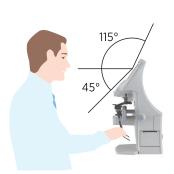



Models and Marking Ink

	White Ink	Red Ink
Cartridge Type	standard	optional
Holder Type	_	standard

^{*} For availability in your region and other details please refer to the local TOPCON authorized representative.

Contact Lens Measurement

Contact lenses can be easily measured using the CL-300 lens support and ring accessories that come standard with the instrument.

Fast, Easy to Load Printer

Loading the thermal paper into the printer has never been easier. Just drop the printer paper into the receptacle and close the cover. The quick and quiet printer is ready to go. This improvement saves time and enhances workflow.



Printout Sample

TOPCON CL-300				
$\langle R \rangle$	S	С	Α	
	+0.00	+0.00	180	
ADD	0.75			
PSM	0.00	OUT	0.00	
〈L 〉	s	С	Α	
	+0.75	+0.25	85	
ADD	0.75			
PSM	0.00	OUT 3.5	50 DWN	
TOPCON Optical Store				

Smart & Smooth Operation

Easy and Fast Operation

The large, newly designed color LCD monitor with a 5.7-inch touch panel simplifies the interpretation of readings even with progressive and multi-focal lenses. The dynamic color measuring graphics follow the measuring process and present each step of the operation with different color, providing better control and ease of use.



In order to do marking, place down the lens retainer to secure the lens and marking is stamped accurately with the equipped

marking cartridge.

Specifications

Measurement Range		Sphere	0 to ±25 D (0.01/0.12/0.25D steps)	
		Cylinder	0 to ±10 D (0.01/0.12/0.25D steps)	
		Axis	0 to 180°(1°step)	
		Addition	0 to + 10D (0.01/0.12/0.25D steps)	
		Prism	0 to 10△	
Cylinder Mode			MIX /-/+	
Prism Mode			X-Y (Orthogonal Coordinates), P-B (Polar Coordinates),No Display	
UV Transmittance			0 ~ 100% (1% / 5% steps)	
Lens Detection			AUTO for single & progressive lenses	
Contact Lens Measurement		ment	Hard / Soft contact lens	
Light Source Refractive Measuremen		ve Measurement	Green LED 540±5nm	
ι	JV Meas	urement	UV LED 365±5nm	
Wavelength Cri	iterion		e-line/d-line with ABBE compensation	
Lens Diameter			φ 5 to 100 mm	
Screen			5.7 inch color touch screen	
PD Measurement Range*		e*	50-90mm	
External Interface			RS-232C , LAN *	
Power Supply			100-240V / 50/60Hz 1.2A (AC Adapter) : AUTO power save	
Dimensions			197(W) × 220(D) ×404(H) mm	
Weight			3.8Kg	

*CL-300PDL model only

Models	CL-300	CL-300PDL
PD measurement	_	included
Printer	included	included
LAN		included

System Chart



Contact Topcon subsidiaries or dealers for system configurations.

*Subject to change in design and/or specifications without advanced notice.

IMPORTANT In order to obtain the best results with this instrument, please be sure to review all user instructions prior to operation.







TOPCON CORPORATION

FC9711081 TOPCON CORPORATION



TOPCON CORPORATION

75-1 Hasunuma-cho, Itabashi-ku, Tokyo 174-8580, Japan. Phone:3-3558-2523/2522 Fax:3-3965-6898 www.topcon.co.jp

TOPCON MEDICAL SYSTEMS, INC.

TOPCON MEDICAL LASER SYSTEMS, INC. Coronado Drive Santa Clara, California 95054 USA ::+1-408-235-8200 Phone:(USA only)+1-888-760-8657 I-408-235-8259 Email: tmlsinfo@topcon.com

TOPCON CANADA INC.
110 Provencher Avenue, Boisbriand, QC J7G 1N1 CANADA
Phone:+1-450-430-7771 Fax:+1-450-430-6457 www.topcor

TOPCON EUROPE MEDICAL B.V.

Essebaan 11; 2908 LJ Capelle a/d Ussel; P.O. Box 2900 AC Capelle a/d Ussel; THE NETHERLANDS Phone: 431 (4)10-4585077 Fax: 431 (4)10-4585045 E-mail: medical@topcon.nl; www.topcon.eu

E-mail: Industria 60; 20037 Paderno Dugnano; (Milano), ITALY
Phone:+39-02-9186671 Fax.+39-02-91081091 E-mail: topconitaly@liscali.it; www.topcon.it

Phone.:439-02-9186671 Fax:-439-02-91081091 E-mail: topconitaly@tiscali.it; ww. DANMARK OFFICE Praestemarkswej 25, 4000 Roskilde, DANMARK Phone:-45-46-327500 Fax:-45-46-327555 E-mail: topcond/epopcondanmark,dk www.topcondanmark.dk IRELAND OFFICE Unit 276, Blanchardstown; Corporate Park 2 Ballycoolin Dublin 15, IRELAND Phone:-353-1897590 Fax:-4353-18293915 E-mail: medical@topcon.ie; www.lopcon.ie

TOPCON S.A.R.L.

TOPCON DEUTSCHLAND G.m.b.H.

Hanns-Martin-Schleyer Strasse 41; D-47877 Willich, GERMANY
Phone:+49-(0)2154-8850 Fax:+49-(0)2154-885177 E-mail:med@topcon.de; www.topcon.de

TOPCON SCANDINAVIA A.B.
Neongatan 2; P.O.Box 25; 43151 Molndal, SWEDEN
Phone:46F.0013-7109240 E-mail:medical@topcon.se; www.lopcon.se

TOPCON ESPAÑA S.A.

HEAD OFFICE
Frederic Mompou 4 Esc. A Bajos 3, 08960 Sant Just Desvern Barcelona, SPAIN
Phone:+34-93-4734057 Fax:+34-93-4733932 E-mail: medica@topcon.es; www.topcon.es

TOPCON (GREAT BRITAIN) LTD.

Topcon House, Kennet Side, Bone Lane, Newbury, Berkshire RG14 5PX United Kingdom Phone: 444(0)1635-551120 Fax: 444(0)1635-551170 E-mail:medical@topcon.co.uk; www.lopco

TOPCON POLSKA Sp. z. o. o. ul. Warszawska 23; 42-470 Siewierz, POLAND Phone:+48-(0)32-6705045 Fax:+48-(0)32-6713405 www.topcon-polska.pl

TOPCON SINGAPORE MEDICAL PTE. LTD.

1 JALAN KILANG TIMOR #09-01 PACIFIC TECH CENTRE SINGAPORE 159303 Phone:+65-68720606 Fax:+65-67736150 E-mail:medical_sales@topcon.com.sg_www.topc

TOPCON INSTRUMENTS (MALAYSIA) SDN.BHD.
No, D1, (Ground Floor), Jalan Excella 2, Off Jalan Ampang Putra,
Taman Ampang Hifir, 55100 Kuala Lumpur, MALAYSIA
Phone:+80-(0)3-42709866 Fax;+80-(0)(3-42709766

TOPCON INSTRUMENTS (THAILAND) CO.,LTD.

77/162 Sinnsathorn Tower, 37th Floor, Krungthonbur Klongsarn, Bangkok 10600, THAILAND Phone: +66(0)2-440-1152-7 Fax: +66-(0)2-440-1158

TOPCON CORPORATION BEIJING OFFICE Block No.9, Kangding Street Beijing Economic-Technological Development Area, Beijing, 100176, CHINA Phone::46-(0)10-6780-2799 Fax:+86-(0)10-6780-2790

TOPCON CORPORATION SHANGHAI OFFICE 14L Huamin Empire Plaza, No.726, Yan-an Xi Road, Shanghai, 200050, CHINA Phone:-86-(0)21-5238-7722 Fax:+86-(0)21-5237-0761 TOPCON CORPORATION BEIRUT OFFICE

TOPCON CORPORATION DUBAI OFFICE P.O.Box 293705,Dubai Airport Free Zone L.I.U. J-12, Dubai, U.A.E Phone:+971-4-299-5900 Fax:+971-4-299-5901