

**The more
Economy Comfort
The better**



OTPB-160



S p a c i o u s . S t r o n g . S t a b l e . S e n s i b l e .

OTPB-160

**The more
Economy Comfort
The better**




The Plantech Motorised Operation Table manufactured and developed exclusively for the Ophthalmic Surgeon is very stable and versatile, It's ergonomic design ensures both surgeon and patient comfort.

The Plantech OPTB 160 is a new generation Motorized Operation Table with new feature and new design. OPTB -160 has a ideal maximum and minimum height ensure Ophthalmologist to perform from convenient, comfortable standing or sitting working position. Its designed allow to clean it overall easily.

Choice of head section give easy access to the patient while allowing maximum leg room and space for other equipment & its paddles. Newly designed main base having castors for free movement and locks & studs are providing stability of patient & table while surgery.

The adjustable cushioned headrest ensures patient comfort and restraint during operation. Specious, Antistatic with deep antistatic polyurethane cushion offers patients excellent comfort and are easy to clean & maintain. The adjustable wrist support provides comfort and stability while surgery. Optimally Positioned head rest with manually controlled up & down motion allows the surgeon easy access to the patient's head/eye position.

Trendelenberg (head low position) and reverse trendelenberg facility ensure safety during emergency.

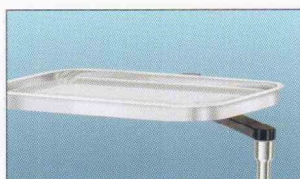

OPERATION TABLE OPTB -160 SPECIFICATION

Maximum Height	: 890 mm
Minimum Height	: 690 mm
Vertical Adjustment Range	: 200 mm
Operational Control	: Foot Switch
Total Length	: 1880 mm
Total Width	: 750 mm
Base Size	: 1000 mm x 460 mm
Base Movement	: Free Movement on Lockable Castor Wheel
Overall Dimension (BodyOnly)	: 1880 mm x 750 mm
Upholstery Length (W/O Headrest)	: 1595 mm
Upholstery Width	: 650 mm
Trendelenberg (Head-Low)	: 29 Degree
Reverse Trendelenberg	: 16 Degree
Auxiliary Out-Put	: 3 No. of 220V AC Power Out Put
Current	: 230 V A/C – 50 Hz – 5 AMP
Fuse	: 5 AMP - 20mm Slow Blow
Motor	: 24V DC Actuator
Weight Capacity	: 200 Kgs
Gross Weight	: 120 Kgs

Additional Accessories

■ Patient Drape Stand ■ I.V. Poles ■ Surgical Instrument Tray ■ Safety Belt

[Since Plantech continually upgrades its technology and know-how, Design & Specifications might undergo positive changes from time to time. The company reserves the right to do this.]